



Added Safety
Unique Design
Adaptable



The Lexa Pro is a unique patient floor lift that features a mast lift system that allows the yoke to move in a vertical or linear direction. The distance between the yoke and mast are constant making it easier to avoid obstacles like bed rails and to allow for easier positioning of patients in a bed or chair. Patients can be lifted from a seated or prone position, even from the floor.



Unique Safety Features

The Lexa Pro's unique safety features set it above all other floor lifts with a braking mechanism incorporated in the lifting arms to avoid a patient free fall. Its unique design not only ensures sway-free lifting but that a sling does not "slide" under a patient when in use. The lifting arm moves seamlessly up and down to lift mast with the use of a hand held control and the tilt yoke can swivel 360° to help with more difficult lifting situations.

allow for a wide variety of slings, including ones to accommodate children, the Lexa Pro is the ideal multi-application floor lift for many environments. The Lexa Pro is equipped with an IBS system (intelligent battery system) that ensures years of trouble-free use.

Quality

Our products are designed and manufactured to meet all the requirements of relevant ISO, UL, IEC and EN standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.



Adaptable

The Lexa Pro's base is spreadable into three positions to accommodate all seating types. Its design provides stability during the lifting situation. With ergonomic push handles, an optional scale and connection options that

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone : +61 1300 132 736
Fax : + 61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories



Slings

Assorted models with clips to accommodate most lifting requirements.




Scale

Must be ordered with lift. Digital reading. Weight displayed in kg or lb. Max. Capacity 350kg/770lb. Part no. L103-1350



Hanger Bar

Available with clip attachment or loop attachment. (Part no. 90337 shown)

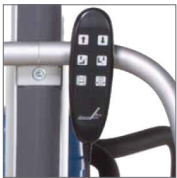

Hand Control

Six functions. Included with floor lift. Replacement Part no. R500-482

Lifting Situations



Use to/from a bed

Use to raise/lower a patient from/to a bed.



Use to/from a chair

Use to raise/lower a patient from/to a chair.



Use to/from the toilet

Use to raise/lower a patient from/to a toilet.

Features




Control Panel

Knee support, mast height (raise/lower) and leg adjustments. Battery low LED indicator.




Dual Wheel Castors

Featuring wheel locking brakes. Easy-swivel design.

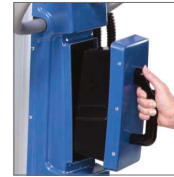

Battery Pack

Easy removal for charging. IBS system.

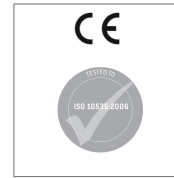


Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Lexa Pro			
Part No.	92501	Lowering Speed	1.8cm/sec 0.7"/sec
Width	67cm / 26.4"	Battery	IBS power pack, 24V/2.9Ah Replacement Part no. R103-1014
Length Length (Lift Arm)	123.9cm / 48.8" 70.2cm / 27.6"	Wheels	Front: 10cm/4" / Optional 7.4cm/3" Rear: 12.5cm/5" / Optional 10cm/4"
Height	Arm Lowered: 175.7cm / 69.2" Arm Raised ¹ : 221.7cm / 87.3"	Hanger Bars (Loop)	2-Point - Part no. 90338 4-Point - Part no. 90339
Height from Floor	13.2cm / 5.2"	Hanger Bars (Clip)	4-Point - Part no. 90337 / 92501
Adjustable Leg (Min - Max)	Inside: 59cm - 120.3cm / 23.2" - 47.4" Outside: 67cm - 133.3cm / 26.4" - 52.5"	Hanger Bars with Scale	4-Point Loop Part no. 92629 4-Point Clip Part no. 90340
Lowering Speed	1.8cm/sec 0.7"/sec	Emergency Lowering	Manual and Electric
Turning Radius	131cm / 51.6"	Weight	70kg / 154.3lb
Sound Level	< 44dB	Weight Capacity	275kg / 606lb

¹ Measurement based on arm mast at highest setting

This product is manufactured in accordance with all relevant ISO standards.

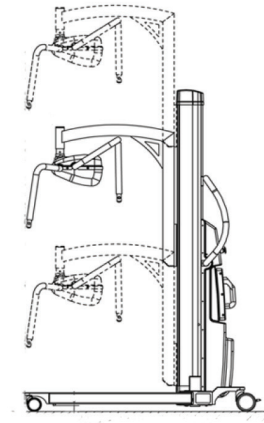

Learn More About This Product



Follow Us on Social Media



facebook.com/humancareusa



instagram.com/humancaregroup



linkedin.com/company/humancareusa



pinterest.com/humancaregroup



youtube.com/channel/UC7fyjNH3zTjRx0P1Df2Da0